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

CE

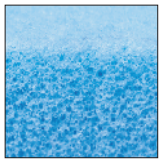WARRANTY 2 LET

VYSOCE PRUŽNÁ PUR PĚNA
- Its mechanical properties (elasticity, flexibility) and high density (40kg/m3) enable positions to be maintained over time, enabling effective treatment.

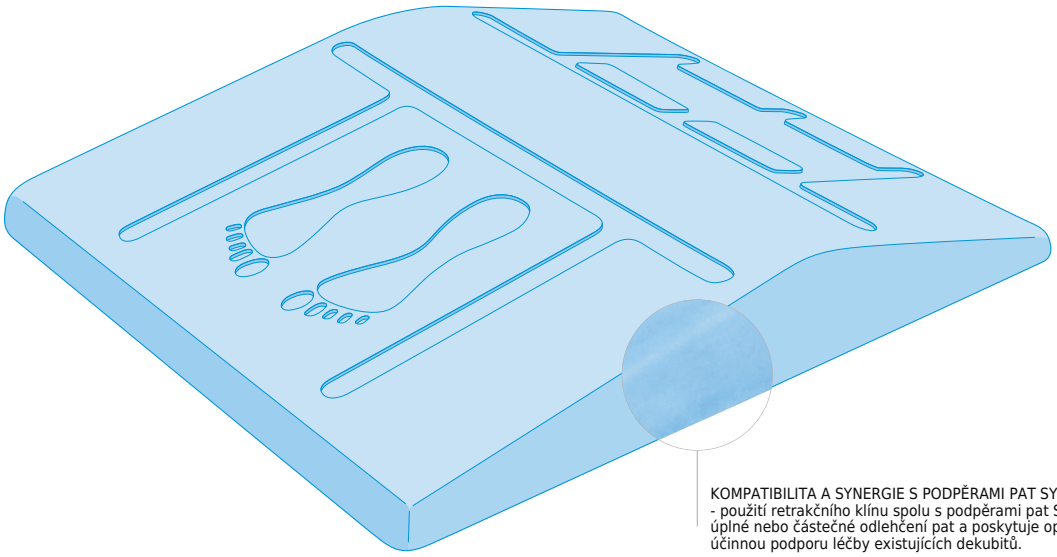EFFICIENCY OF THE SEMI-FOWLEROVA POLOHA
- Facilitates the relief of pressure on the pata and sakrum areas, by transferring pressure to the calves and thighs (zones at low risk of dekubit formation).
- Ensures patient stability and reduces "forward sliding", limiting the effects of tření and střižné síly which play a major role in the process of dekubit formation.

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LEARN MORE

RELEVANCE OF THE SEMI-FOWLEROVA POLOHA

Stiffness and oppositional hypertonia in geriatrics, factors to be considered in postural strategy for dekubit prevention

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Less 5° 5 to 10° More 10°

Angle measurements for articular relaxation of the knees in the dorzální dekubitus position 56 measurements, average age of 82 (fig. 1)

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SYSt’aM® P910L / 30° POLOHOVACÍ KLÍNY

POSITIONING KLÍN MADE OF MEMORY FOAM AND A STABILISATION NUB MADE OF VYSOCE PRUŽNÁ PUR PĚNA

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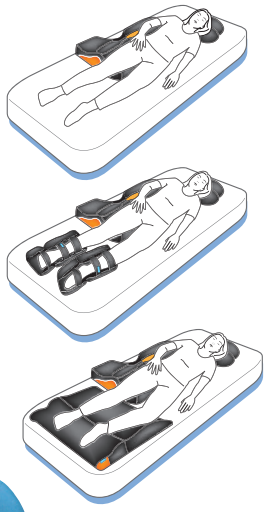
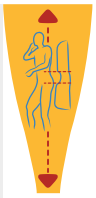
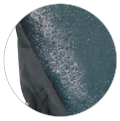
ORIENTAČNÍ ŠTÍTEK
- Enables perfect positioning of the sakrum pressure-release system on the right physical area of the patient: controlled and increased effectiveness of use.
- Facilitates positioning for care staff.

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
CARRYING HANDLE
- Placed on the lower face.
- Facilitates handling and transfer from one room to another.

OZNAČENÍ	ROZMĚR D x Š x V (cm / inches) KÓD POLOŽKY	
P910L / 30°- Positioning wedge	P910L1HW	99 x 55 x 15 cm / 39 x 21,6 x 6’’
P910L*D* / Half 30°- Positioning wedge	P910LD1HW	99 x 28 x 15 cm / 39 x 11,2 x 6’’
P910L*G* / Half 30°- Positioning wedge	P910LG1HW	99 x 28 x 15 cm / 39 x 11,2 x 6’’


CE WARRANTY 2 LET

05/2019

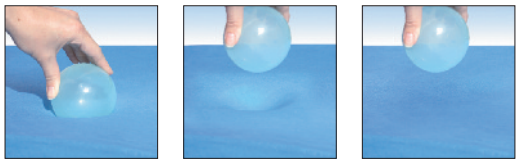TVAROVANÁ VISKOELASTICKÁ PAMĚŤOVÁ PĚNA
- Precise moulding of the body and increase of the body surface in contact s the mattress;
- Better pressures distribution: reduction of transcutaneous pressures on areas at high risk and facilitated blood circulation;
- Improved comfort and stability of the patient;
- Skin effect obtained through the moulding process: protection of the foam against external aggressions (tear, crumbling);
- With a very high density to prevent deformation and sagging effects.



OPORA DOLNÍ KONČETINY V OBLASTI STEHNA
- Contributes to relieving pressure on the sakrum area by transferring opora to a lower riziková oblast.
- Enhances patient comfort by avoiding rotation of the vertebral column, contributing to more effective use.
- Consequently lowers pressure exerted on the pelvic area.
- Prevents contact between the knees (a painful area posing the risk of dekubit s prolonged contact).

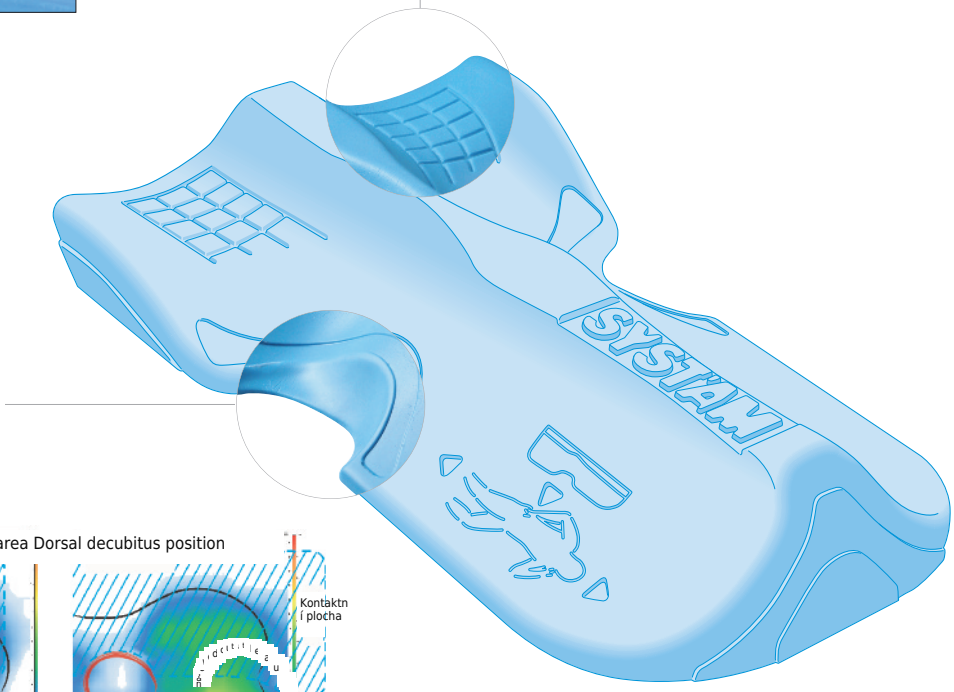VYBRÁNÍ V OBLASTI KŘÍŽOVÉ KOSTI
- Relieves pressure on the riziková oblast.
- Contributes to reduce painful contacts.
- Increases the proven efficiency of the position at 30°.

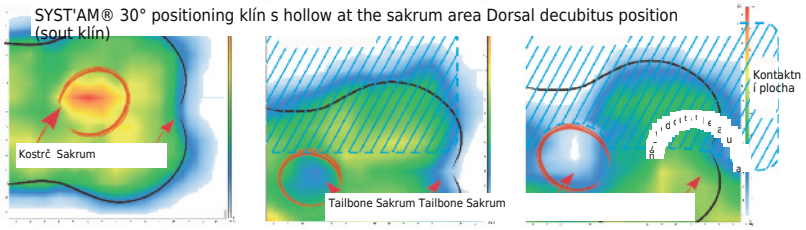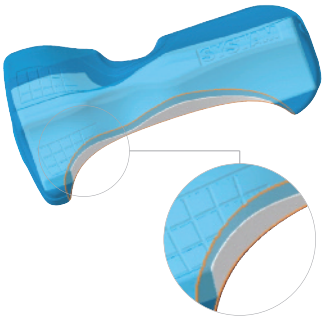SYST'AM® 30° positioning klín s hollow at the sakrum area Dorsal decubitus position (sout klín)


STABILIZAČNÍ VLOŽKA The SYST'AM® 30° klín has been designed using multi-material stratification combining a central anatomical nub foam moulded to an very soft viscoelastic surface.
- The nub is a key feature for the effective use of the device, enabling the 30° position to be maintained over time by acting as a opora structure.
- The nub enables to have the upper part made of exceptionally soft viscoelastic foam, ensuring comfort and a consistent reduction in pressure.
- The ergonomic shape of the insert respects the body curves and fade out accordingly at areas at risk.


30° SYMETRICKÝ KLÍN
- Can be used on both the right and left side.
- Facilitates handling during position reversal.

30° ANATOMICKÝ KLÍN
- Respects and moulds to bodily curves.
- Contributes to patient comfort and effective use.

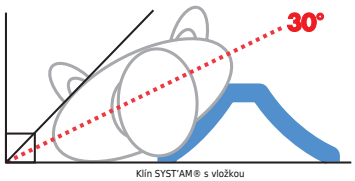ZMENŠENÁ ŠÍŘKA KLÍNU
- Allows the patient to be perfectly centred on the bed.
- Aims to avoid accidental contact of the patient s the barriers of the bed.
- Will not become lodged between the barriers of the bed.

Klín SYST'AM® s vložkou

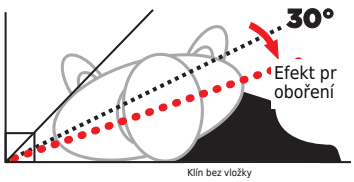Klín bez vložky

LEARN MORE

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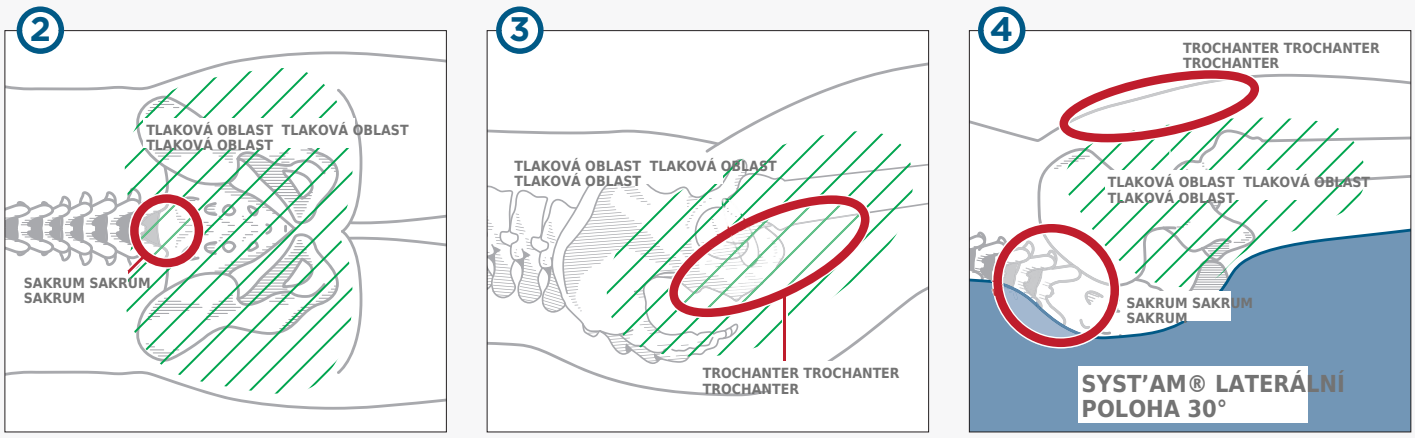
In general, dekubit prevention and treatment help is fundamentally based on the principle of reducing pressure on tissue, by increasing the area of surface contact or transferring opora to low-rizikové oblasti.
- Complementary to this, the majority of medical departments include changes of position as part of their strategy of dekubit prevention for patients at high or very high risk, s a lateralisation of 90° in practice.
- While this position effectively removes pressure on the sakrum region, it creates significant additional pressure on the greater trochanter area.

The use of such a position is often responsible for the development of multiple dekubity, worsening the vital prognosis.
- Studies show that patients at risk of developing dekubity do not regain full oxygenation of tissue between different phases of postioning.
- Frequent changes in position are therefore not sufficient in preventing the risk of trochanterian dekubity in the lateral position at 90°.

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
The semi-lateral position at 30° is widely described and recommended by several authors. The French national agency of accreditation (ANAES) advocates the use of the lateral-oblique decubitus position at 30° to replace the lateral decubitus position at 90°, given the resulting risk of trochanterian dekubits.
- The principle behind lateral positioning at 30° is the protection of rizikové oblasti (essentially the sakrum and trochanters) by shifting pressure to zones at low risk of dekubits, sout bony protuberance and well vascularised (the postero- external side of the pelvis).
- An analysis of the effect of position on distribution of pressure conducted by Delfloor et al. concludes that the semi-lateral decubitus position at 30° shows distinctly lower maximum pressure than otherweise recorded s any other position.
- Multiple studies, particularly Seiler et al, have been conducted to compare the condition of tissue oxygenation at the sakrum and the trochanters in relation to position. Measurements were recorded in the dorzální dekubitus position and in positions of 90° and 30°.
- The results confirm the appearance of full or severe skin hypoxia at the sakrum in the dorzální dekubitus position, and also showed trochanterian skin hypoxia at 90°.

In the semi-lateral decubitus position at 30°, sacro-gluteal TcPo2 remains close to the value at rest, regardless of which opora is used.
- The studies also show the preservation of skin oxygenation at the trochanter in the semi-lateral decubitus position at 30°. Conclusion:
- Positioning at 30° consequently enables the protection of body parts which are commonly positioned in the dorzální dekubitus position (essentially the sakrum) from the risk of dekubity, by avoiding the transfer of this risk to the trochanters (positioned at 90°).
- At the same time, the studies show that the benefits of positioning at 30° are clearly not dependent on the nature of dekubit prevention mattress used.

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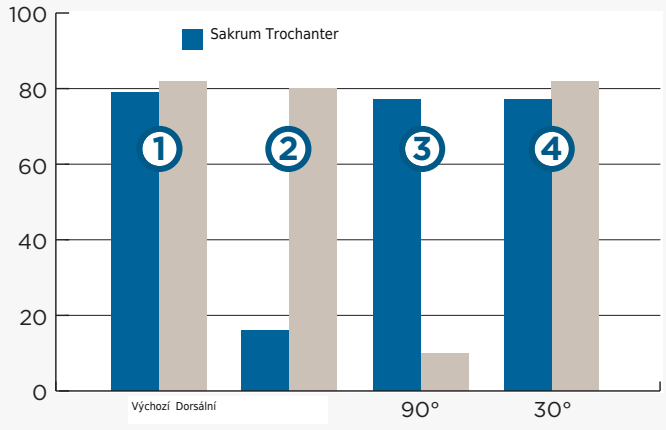Lateral positioning at 30° is difficult to sustain. Seiler describes the semi-lateral position as the maintenance of a position s the back angled at 30°, s the leg of the raised side semi-flexed at the hip and the knee (an essential condition for the sakrum to be completely relieved of pressure).
- In practice, this position is difficult to achieve, especially when the patient is awake and has limited mobility.
- A number of bed oporas are currently used to perform patient lateralisation (polštářs, pillows, bolsters, blankets, foam shapes etc...). However, such oporas are often shaped sout anatomical curves and consistency is rarely satisfactory, therefore they consequently do not enable the 30° position to be correctly controlled over time.
- As a result, ineffective use and a decline in the angle were observed s these types of opora over time.
- In addition, the lack of opora on the raised side of the thigh led to the rotation of the vertebral column, causing great discomfort to the patient and consequently exerting great pressure on the pelvic area.

VLIV NA OBLAST PATY

A literature review of the impact of the semi-lateral position at 30° on pata pressure showed in all cases, that this area suffered less pressure than in the decubitus dorsal position.
- Despite the opora of the thigh limiting pressure exerted on the pata on the raised side in the semilateral decubitus position, the addition of a pata opora can prove complimentary as a strategy to facilitate high levels of prevention.

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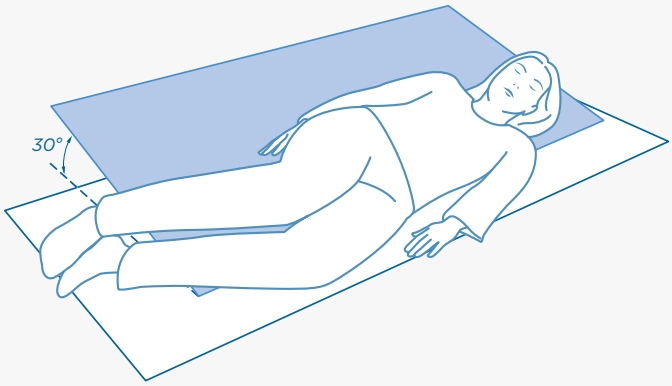
Specially designed for lateral positioning at 30° (in contrast to pillows, bolsters, foam blocks etc...).

Anatomical shape in the form of physical curves (in contrast to pillows, bolsters, foam blocks etc...).

Designed to facilitate total relief of pressure on the sakrum
- by oporaing the thigh of the raised side (patient comfort, reduction in pressure on the pelvic area)
- by a hollowed out area of the device located at the sakrum (sakrum completely relieved of pressure).

Ability to maintain the 30° position over time, coupled s a comfortable surface and high capacity for pressure reduction (the benefit of an insert contributing to stabilisation and a surface made of viscoelastic memory foam).

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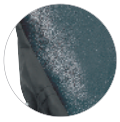
P912L is entirely protected by a potah and has a particularly comfortable attachment system that is easy to adjust:
- the attachment system offers sustainable opora maintaining the device in the right position,
- the system is easy to install.

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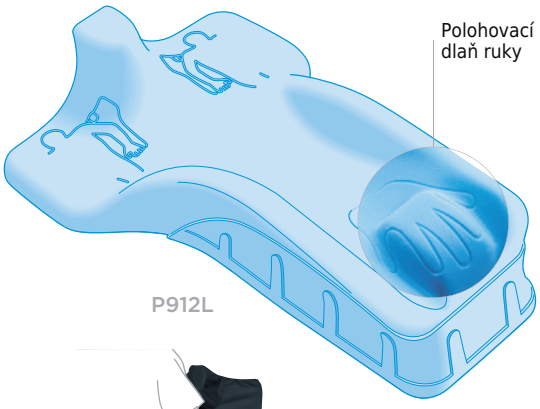
ANATOMICKY TVAROVANÝ KLÍN
- The anatomical shape respects the natural curve of the arm in order to guarantee patient comfort and maintain the position of the horní končetiny over time, for more effective use.
- Specially designed shape built in accordance s recommendations for horní končetiny declining positions.
- Device is easy for the medical team to use.

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V POLOZE VLEŽE
- The extension of the elbow avoids restricting vein perfusion in the forearm.
- Allows access for a drip to be put in place.

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
COVER WITH ATTACHMENT SYSTEM Comfortable: Unlike most abdukce systems using belts and Velcro fasteners, which are uncomfortable and cause irritation.
- Easy to adjust: is easy to install for caring staff.
- Offers sustainable opora maintaining the legs in the right position.
- Unlike most systems attached to dolní končetiny using belts and Velcro fasteners, which are uncomfortable and cause irritation, the SYST'AM® P911L is entirely protected by a potah and has a particularly comfortable attachment system that is easy to adjust:
- the attachment system offers sustainable opora maintaining the device in the right position,
- the system is easy to install.

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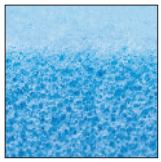
CE WARRANTY 2 LET

VYSOCE PRUŽNÁ PUR PĚNA
- No sharp edges, thus reducing the zones subject to střižné síly and the resulting discomfort.
- High-resiliency foam for optimal comfort in pressure areas.

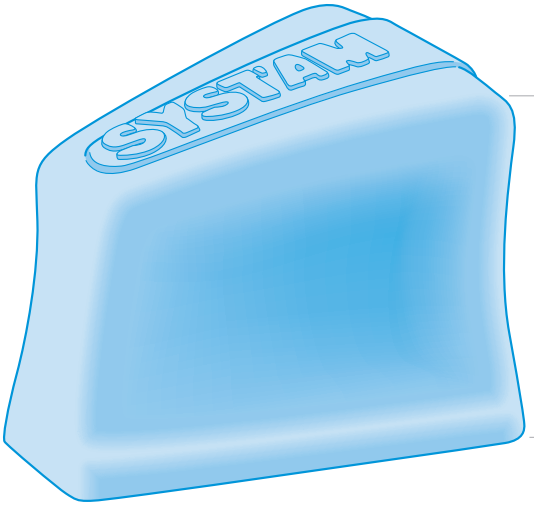ANATOMICALLY SHAPED ABDUKČNÍ KLÍN
- Pre-shaped sides respecting the morphology of the thighs.
- The shape fits the natural angle between the dolní končetiny (distal part finer than the proximal part) thus ensuring optimal distribution of the pressure zones.
- Fine curved tip to help carers to position the klín effortlessly.

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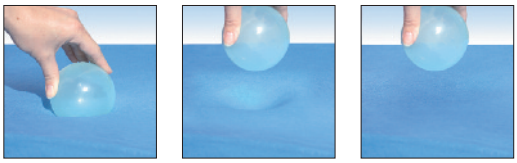
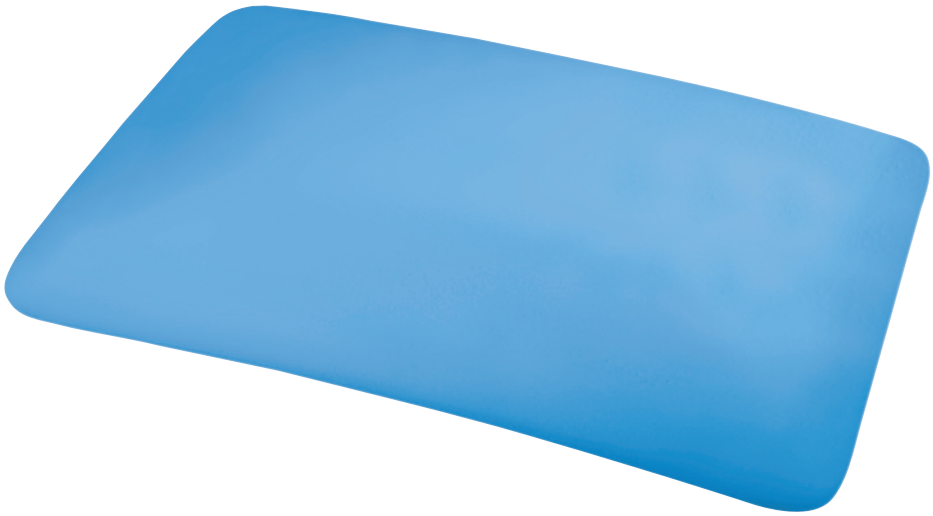
Čištění potahu POLYMAILLE®:

INDIKACE

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

VISKOELASTICKÁ PAMĚŤOVÁ PĚNA
- Defined body moulding and increased body surface area in contact s the foam.
- More effective distribution of pressure: reduction of transcutaneous pressure on highrizikové oblasti, facilitated blood flow.
- Increased comfort and more effective maintenance of patient position.

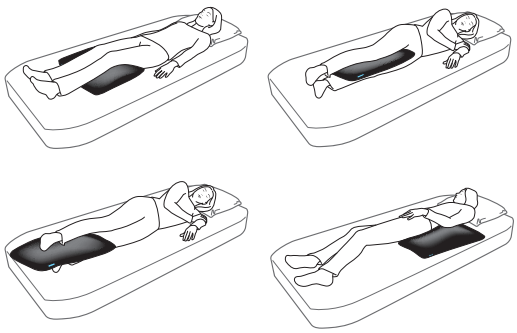This is the perfect oporaing aid: it helps opora the patient, relieves the pressure on the patas, protects the internal lateral knee areas, ankles...
- This polštář can also be used as postural help in Rest chair, wpatachair and shell seat.
- Facilitates positioning for care staff.
- Used by itself or s other devices.

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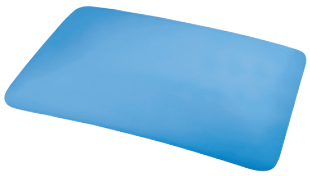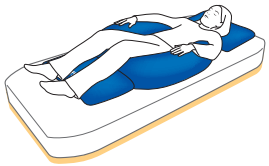
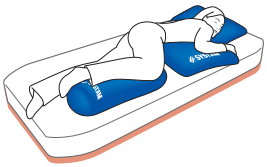
This is the perfect oporaing aid: it helps opora the patient, relieves the pressure on the patas, protects the internal lateral knee areas, ankles and can be used as pillow.
- Used by itself or s other devices.
- Comes in two sizes, Small and Standard (see sizes on last page).



ABDUCTION BLOCK FOR KNEES (Kód položky P9703B1HW / P9753F1HW)
- The simultaneous use of P9702 and P9703 polštářs achieves a comfortable semi- Fowler position for the patient (back and pelvis at 30°) as well as a dorsal position by firmly holding the pelvis and abdukované areas in place thanks to the abdukce semi-circle polštář P9703.
- Reduces pressure on the internal malleolus and the inner sides of the knees.

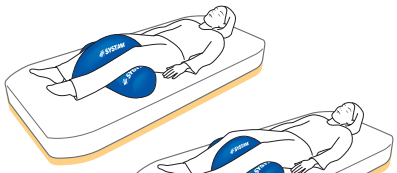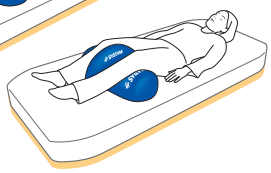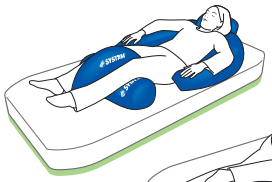
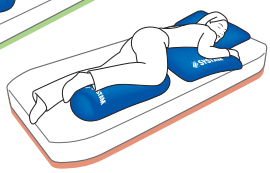
CYLINDRICAL POLŠTÁŘ (Kód položky P9702B1HW / P9752F1HW)
- Can be used in various ways, such as placed under the knees to maintain dolní končetiny in a Semi-Fowler position.
- When used s the abdukce block (P9703), it firmly holds abdukované limbs in position and protects the inner sides of the knees.


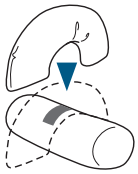
SUCHÝ ZIP
- jednoduchý a účinný.

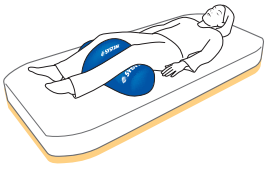RING POLŠTÁŘ (Kód položky P9704B1HW / P9754F1HW)
- Comfortably oporas the head both in sitting or semi-reclined position.
- Easy-to-use, self-gripping strip.

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

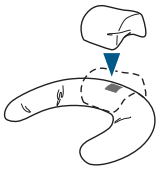
Holds the pelvis and upper body firmly in place when patient in Semi-Fowler or dorsal position. Relieves pressure on elbows. Fitted s a strap for extra opora P9706 enabling to maintain the head in position.
- Can also be used on a wpatachair.
- SYST'AM® exclusivity: equipped s the EasyShape® system that enables to maintain a position that is adapted to the patient's morphology.

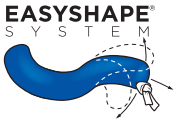P9708B / P9758F
SUCHÝ ZIP
- jednoduchý a účinný.


LATERAL POSITIONING POLŠTÁŘ (Kód položky P9707B1HW / P9757F1HW)
- Designed s integrated pillow for optimal comfort.
- Ensures comfortable positioning of the knees, pelvis and upper body when in lateral position.
- SYST'AM® exclusivity: equipped s the EasyShape® system that enables to maintain a position that is adapted to the patient's morphology.
- Available in one size, compatible s all patient sizes.


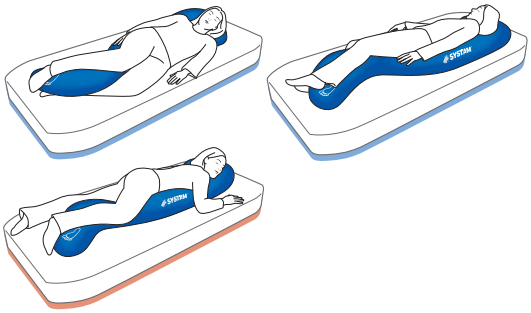
EASYSHAPE®
SYSTEM


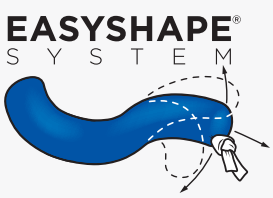
EASYSHAPE®: EXCLUSIVE SYST'AM® ADJUSTMENT SYSTEM
- The straps attached to the half-moon polštář potah (P9705B / P9755F) and to the lateral positioning polštář potah (P9707B / P9757F) provide an easy way to adjust and to maintain the device in a position that is adapted to the patient's morphology, thereby increasing the patient's comfort and enhancing the ergonomics of the opora.
- By pulling on the straps the curvature of the polštář increases and the shape is maintained by tying the straps.
- Simple and easy to use, the EasyShape® system enables to adjust the device to better adapt it to the patient's morphology and to ensure that the position is maintained for an extended period of time.

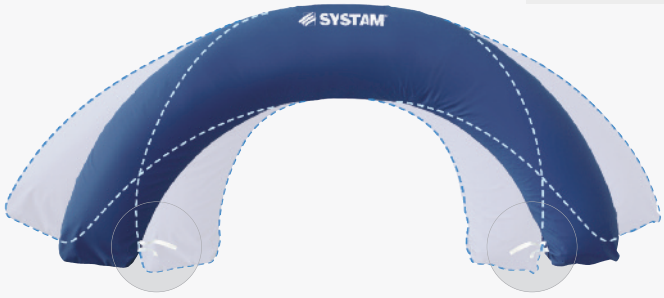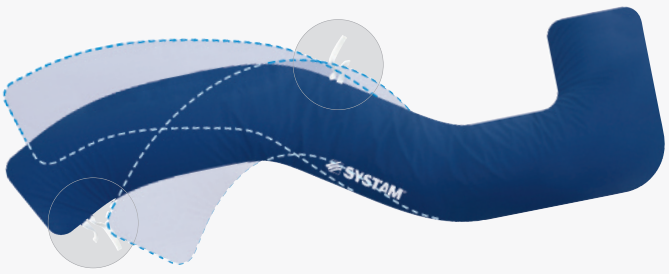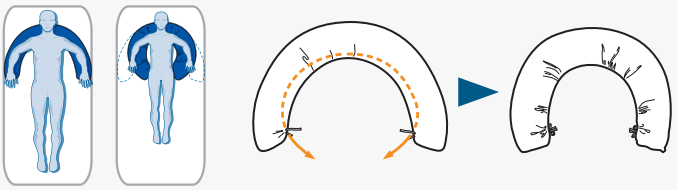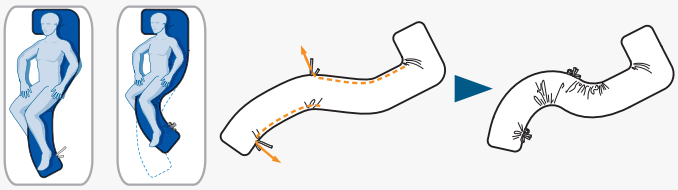
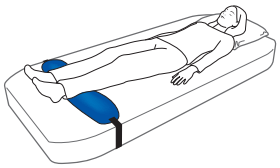
CROSS PODLOŽKA WITH ADJUSTABLE STRAP

- Enable hips and knees positioning as well as full pressure relieve of the patas.
- Adjustable strap to fix on the mattress.

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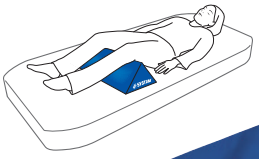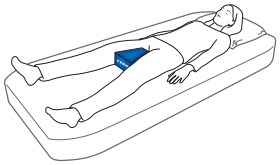
- Pressure reduction on the occipital area.

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Disinfect using cold-application spray products (quaternary ammoniums, glutaraldehydes,etc).

LEARN MORE

BIBLIOGRAFIE

PODPĚRY PATY

[1] ABU-OWN A, SOMMERVILLE K, SCURR JH, COLERIDGE S . Effects of compression and type of bed surface on the microcirculation of the heel. Eur. J. Vasc. Endovasc. Surg. 1995 apr, 9(3): 327-34. [2] ALLEN V, RYAN DW, MURRAY A . Air-fluidized beds and their ability to distribute interface pressures generated between the subject and bed surface. Physiol. Meas. 1993 aug, 14(3):359-64. Barton aa, barton m . Drug-based prevention of pressure-sores. Lancet 1976,ii: 443- 4. [3] BLASZCZYK J, MAJEWSKI M, SATO F. Make a difference: standardize your heel care practice. Ostomy Wound Manage. 1998 may, 44(5): 32-40. [4] BURDETTE – TAYLOR SR, KASS J. Heel ulcers in critical care units: a major pressure problem Crit Care Nurs Q 2002 Aug: 25(2): 41-53. [5] CHENEY AM. Portrait of practice: a successful approach to preventing heel pressure ulcers after surgery. Decubitus . 1993 jul, 6(4): 39-40. [6] CHENEWORTH CC, HAGGLUND KH, VALMASSOI B, BRANNON C . Portrait of practice: healing heel ulcers. Adv. Wound Care. 1994 mar, 7(2): 44-8. [7] DE KEYSER G, DEJAEGER E, DE MEYSTH, EDERS GC. Pressure-reducing effects of heel protectors. Adv. Wound Care. 1994 jul, 7(4): 320-2. [8] FONTAINER R, RISLEY S, CASTELLINO R. A quantitative analysis of pressure and shear in the effectiveness of support surfaces. J ; Wound Ostomy Continence Nurs. 1998 sept, 25(5): 233-9. [9] CRAFF MK, BRYANT J, BEINLICH N . Preventing heel breakdown. Orthop. Nurs. 2002 sept oct, 19(5): 63-9. [10] HUNTER SM, LANGEMO DK, OLSON B, HANSON D et al. The effectiveness of skin care protocols for pressure ulcers. Rehabil. Nurs. 1995, 20: 250-55. [11] JORDAN MM, CLARK MO ; Report on the incidence of pressure sores in the patient community of the greater glasgow health board area, 21st january, 1976. Bioengineering, unit, university of strath-clyde, and the greater glasgow health board, 1977. [12] MAKLEBUST J, SIGGREEN MY, MONDOUX L Pressure relief capabilities of the Sof. Care bed and the Clinitron bed. Ostomy Wound Manage 1988 Winther , 21: 32, 36- 41, 44. [13] MEEHAN M. National pressure ulcer prevalence survey. Adv. Wopund Care. 1994, 7 27-38. [14] PARNHAM A. Interface pressure measurements during ambulance journeys. J. Wound Care. 1999 jun, 8: 279-82. [15] PINSUR MS, SCHUMACHER D, REDDY N. Preventing heel ulcers: a comparison pf prophylactic body-support systems. Arch Phys Rehabil. 1991 jun, 72: 508-10. [16] TOURTUAL DM, RIESENBERG LA, KORUTZ CJ, et all. Predictors of hospital acquired heel pressure ulcers. Ostomy Wound Manage. 1997 oct, 43: 24-8, 30,32-4. [17] TYMEC AC, PIEPER B, VOLLMAN K, A comparison of two pressure-relieving devices on the prevention of heel pressure ulcers. Adv. Wound. Care. 1994 jan feb 7, 10: 39-44.

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[1] Delflor T. The effect of position and mattress on the development of pressure sores. In: Verpleegkunde 1997, 12:140-9. [2] Collin D, Abraham P, Perault L Comparison of 90 degrees and 30 degrees laterally inclined positions in the prevention of pressure ulcers using transcutaneous oxygen and carbon dioxide pressures. In: Advances in Wound Care 1996, 9:35-8. [3] Bader D.L., Hawken M.B. Ischial pressure distribution under the seated person. In: Pressure sores. Clinical practice and scientific approch 1993, 17:223-33. [4] Bouttereux B., Collin D. L'intérêt du décubitus semi-latéral pour décharger l'appui talonnier. In: L'escarre. Revue officielle de l'association PERSE. 2003, 18:15. [5] Sachse R.E., Fink S.A., Klitzman B. Comparison of supine and latéral positioning on various clinically used support surfaces. In: Ann. Plast. Surg. 1998, 41:513-8. [6] Seiler W.O., Allen S., Stahelin HB. A new investigative method using transcutaneous oxygen tension measurement. In: J. of Am. Geriat. Soc. 1983, 31:786-9. [7] Seiler WO., Chapuis A., Stahelin HB., Dolfus P. Prevention of decubitus ulcer. "Supersoft" support and 30 degrees lateral-oblique positioning. Reducing the pressure as much as possible. In: krankenpfl. Soins Infirm. 1984, 77:51-3. [8] Seiler W.O, Allen S.,Stahelin H.BInfluence of the 30 degrees laterally inclined position and the "super soft" 3 pieces mattress on skin oxygen tension on areas of maximum pressure implications for pressure sore prevention. In: Gerontology 1986, 32:158-66. [9] Schubert V., Heraud J. The effect of pressure and shear on skin microcirculation in elderly stroke patients lying. in supine or semi-recumbent positions. In: Age Ageing. 1994, 23:405-10. [10] Seiler W.O., Stähelin H.B. Skin oxygen tension as a function of imposed skin pressure: implication for decubitus ulcer formation. In: J Am Geriatric Soc. 1979, 27:298-301. [11] Newson B.A., Pearcy M.J., Rolfe P. Skin surface PO2 measurement and the effect of externally applied pressure. In: Arch Phys Med Reabil. 1981, 62:390-92.

UPPER MEMBER POSITIONING KLÍN

[1] ANAES. Prise en charge initiale des patients adultes atteints d'accident vasculaire cérébral. [2] Corley FG, Schenck RC. Fractures of the hand. Clin. Plast. Surg. 1996;23(3):447-62. [3] Clark S. Brown DM, Elevation crutch in the treatment of the edematous hand. Am J Occup Ther. 1978 May-Jun;32(5):320-1. [4] Eysette M, Seyez JF, Boisson D, Girard R. l'épaule de l'hémiplégique: aspects cliniques à propos de 100 cas. Lyon Méd. 1984, 252:7-13. [5] Giudice ML. Effects on continuous passive motion and elevation on hand oedema. Am. J. Occup. Ther. 1990;44(10):914-21. [6] Kotzki N, Pelissier J, Dussotoit C et al. Techniques de prévention du syndrome algodystrophique: évaluation d'un protocole d'installation au lit. Ann Réadapt. Méd. Phys. 1991,34:351-5. [7] Luce EA. The acute and subacute management of the burned hand Clin Plast Surg. 2000 Jan;27(1):49-63. [8] MacDermid JC. Hand therapy management of intra-articular fractures with open reduction and pi plate fixation: a therapist's perspective. Miles W. Soft tissue trauma. Hand Clin. 1986 feb;2(1):33-43. [9] Olivera MA, Araujo JF. Carpal tunnel syndrome: retrospective study of 106 cases ans complications. Arq Neuropsiquiatr. 1993;51(4):519-24. [10] Reid DA. Post-operative complications in hand surgery. Hand. 1975 Jun;7(2):115-19. [11] 3Sorenson MK. The edematous hand. Phys. Ther. 1989;69(12):1059-64. [12] Silver L. Burns of the hand and wrist. Current concepts in diagnosis and treatment. Orthop Rev. 1987 Jun;16(6):394-400. [13] Tromel M.F., Leclerc J. Rééducation de la main brulée Ann. Médit. Burns Club 1994,11:3. [14] Uzelac O, Nikolic Z. Postoperative treatment in the surgery of acute hand injuries. Acta. Chir. lugosl. 1982;29(2):323-8. [15] Vasudevan SV, Melvin JL. Upper extremity control: rationale of the techniques. Am. J. Occup. Ther. 1979;33(8):520-3.